

Bepotastine besilate

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK®ID [4.6 x 250mm, 5µm]
Part No: 84325

Mobile phase: n-Hexane : EtOH : 1,4 Dioxan : TFA : DEA
(80/10/10/0.1/0.1, v/v/v/v/v)

Flowrate: 1.0 mL/min


Column Oven Temperature: 25 °C

UV detection: 260 nm

Diluent: Ethanol : Solution A (2:8)

Solution A : n-Hexane : EtOH : 1,4 Dioxane (8:1:1)

Concentration: 2.0 mg/mL.

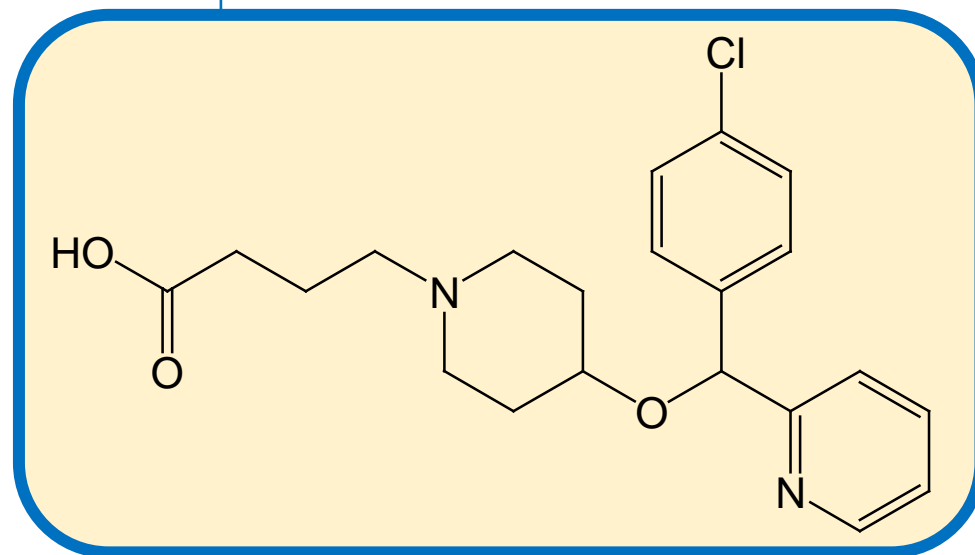

SEPARATION PERFORMANCE

tR ₁ : 8.80	N ₁ : 8491	T ₁ : 1.3	
tR ₂ : 10.74	N ₂ : 3240	T ₂ : 1.9	Rs: 3.42

