

Besifloxacin

CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I(-) [3.0 x 150 mm, 5 µm]

Part No: 54784

Mobile phase: aq. Perchloric acid (pH=1.0):ACN (60/40,v/v)


Flowrate: 0.5 mL/min

Column Oven Temperature: 15 °C

UV Detection: 295 nm

Diluent: Mobile phase

Sample Concentration: 0.4 mg/mL

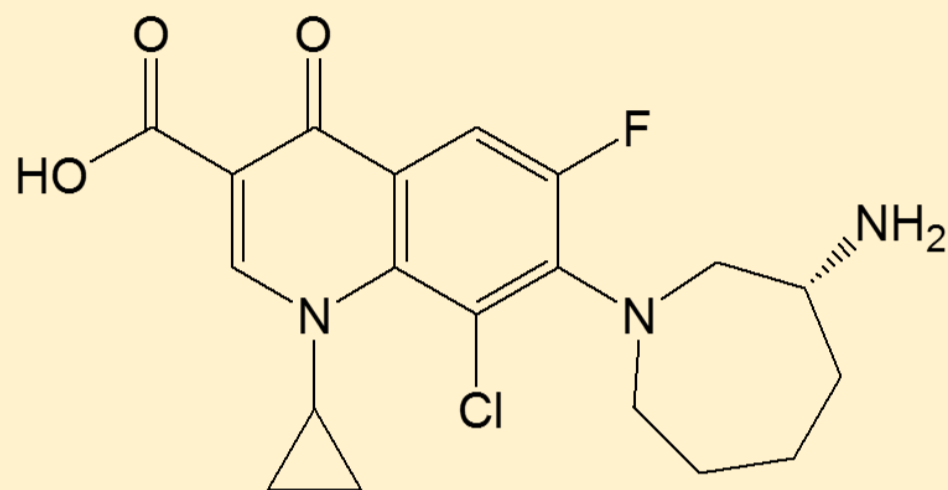

SEPARATION PERFORMANCE

tR ₁ : 19.7	N ₁ : 5344	T ₁ : 1.1	
tR ₂ : 22.3	N ₂ : 4707	T ₂ : 1.4	Rs: 2.25

