

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Brivaracetam

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 150 mm, 5 µm]

Part No: 87324

Mobile phase: aq. 10 mM ammonium bicarbonate : ACN (60/40, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 215 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

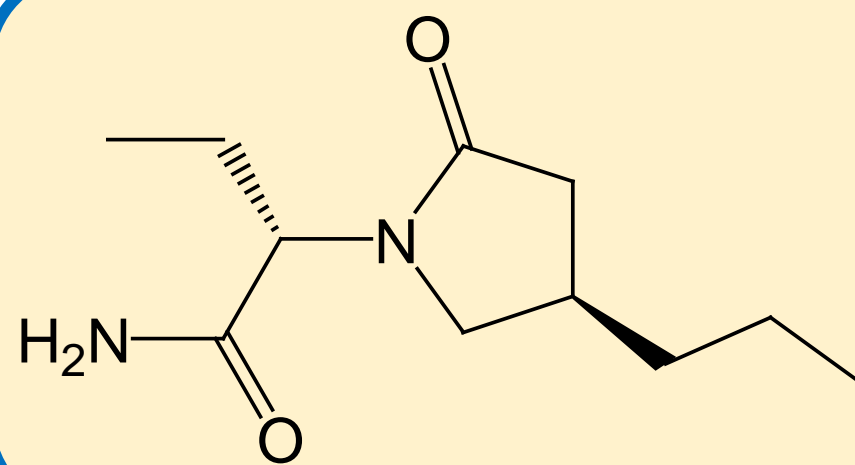

SEPARATION PERFORMANCE

tR ₁ : 5.2	N ₁ : 6700	T ₁ : 1.2	
tR ₂ : 7.4	N ₂ : 6200	T ₂ : 1.2	Rs: 7.1

