

Crizotinib

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ID [4.6 x 250 mm, 5 µm]
Part No: 84325

Mobile phase: n-Hexane : 2-Propanol : Methanol : iso propyl amine
(50/30/20/0.3, v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 265 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

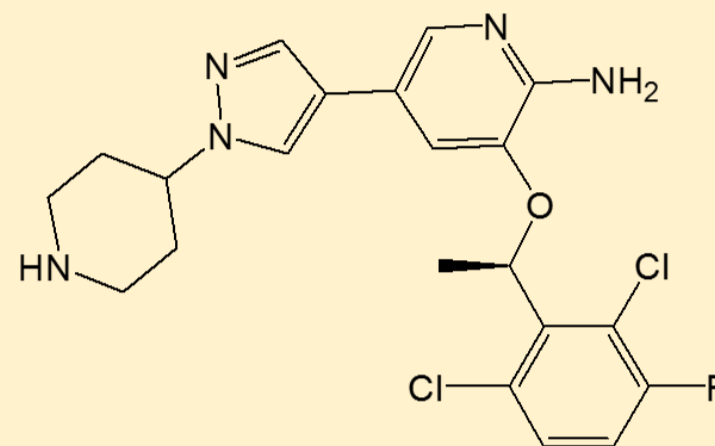

SEPARATION PERFORMANCE

t _{R1} : 11.2	N ₁ : 6247	T ₁ : 1.5	
t _{R2} : 14.4	N ₂ : 6411	T ₂ : 1.4	Rs: 4.9

