

Efavirenz

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IC-3 [4.6 x 150 mm, 3 µm]
Part No: 83524

Mobile phase: 5mM ammonium bicarbonate in
(H₂O : Methanol 20/80, v/v)


Flowrate: 0.5 mL/min

Column Oven Temperature: 25 °C

UV Detection: 246 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

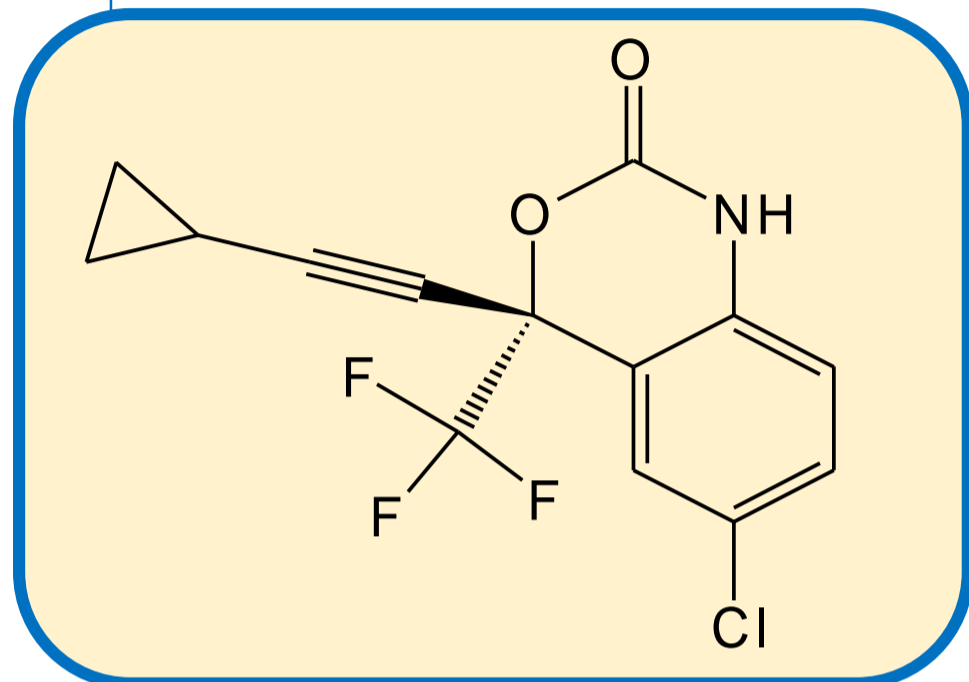

SEPARATION PERFORMANCE

t_{R1} : 6.2	N_1 : 13464	T_1 : 1.3	
t_{R2} : 7.5	N_2 : 12208	T_2 : 1.3	R_s : 5.2

