

Efinaconazole

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK®IG [4.6 x 250 mm, 5 µm]

Part No: 87325

Mobile phase: Methanol (100%)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV detection: 210 nm

Diluent: Methanol

Concentration: 1.0 mg/mL.

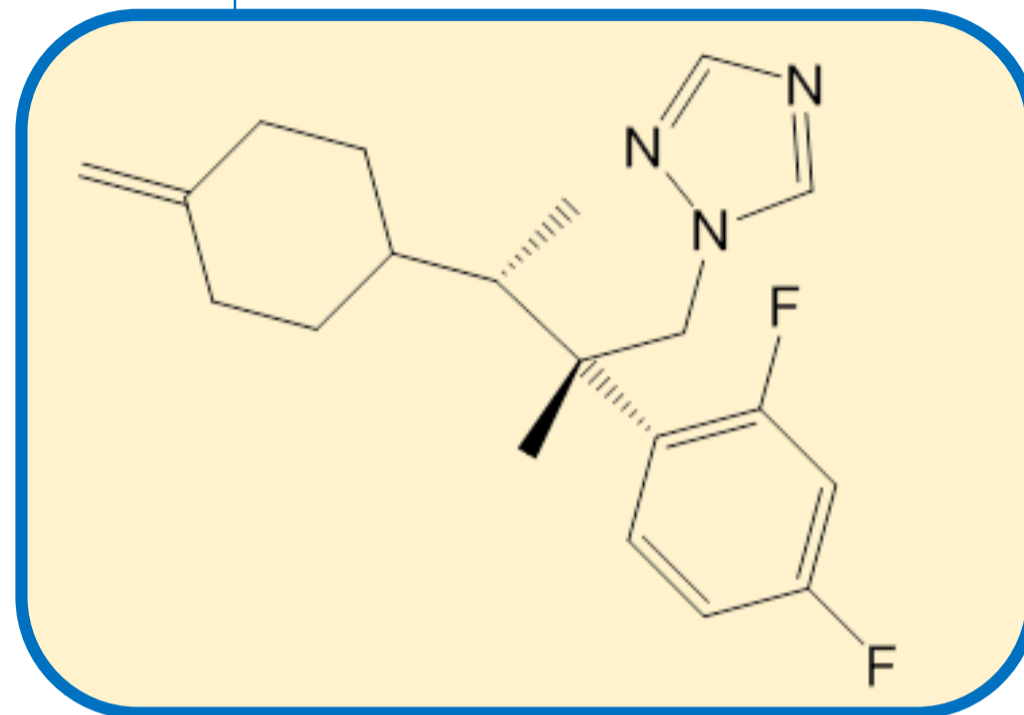

SEPARATION PERFORMANCE

t_{R_1} : 5.45	N_1 : 9124	T_1 : 1.2	
t_{R_2} : 8.68	N_2 : 5826	T_2 : 1.2	R_s : 9.44

