

Eletriptan

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]
Part No: 87525

Mobile phase: n-Hexane : Ethanol : Methanol : DEA
(70/10/20/0.1, v/v/v/v)


Flowrate: 0.8 mL/min

Column Oven Temperature: 20 °C

UV Detection: 225 nm

Diluent: Ethanol

Sample Concentration: 1.0 mg/mL

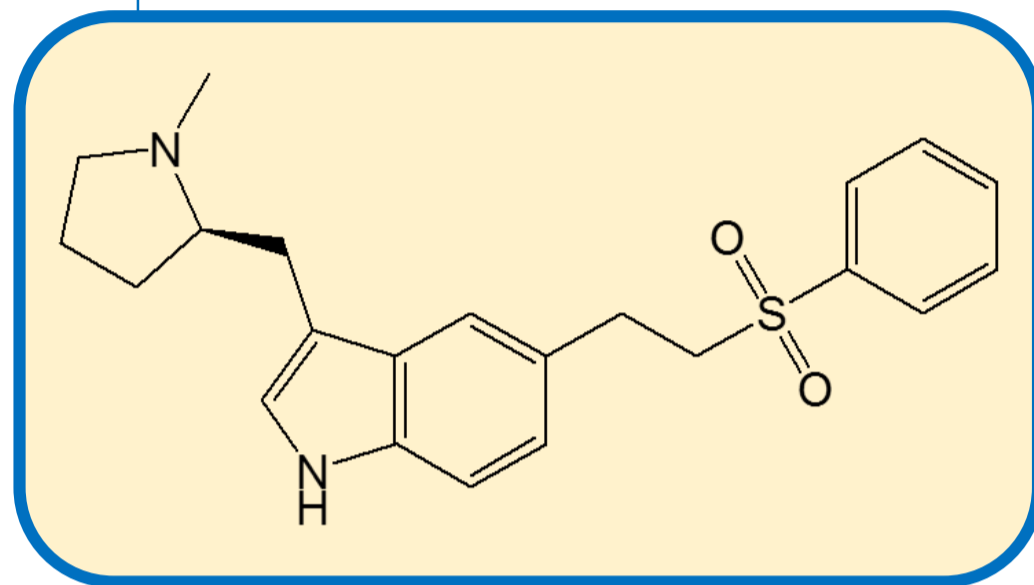

SEPARATION PERFORMANCE

tR ₁ : 13.0	N ₁ : 10181	T ₁ : 1.3	
tR ₂ : 14.1	N ₂ : 8769	T ₂ : 1.4	Rs: 2.0

