

## Levomilnacipran

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: n-Hexane : EtOH : MeOH : EDA (80/10/10/0.1; v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 225 nm

Diluent: Ethanol (100 %)

Sample Concentration: 1.0 mg/mL

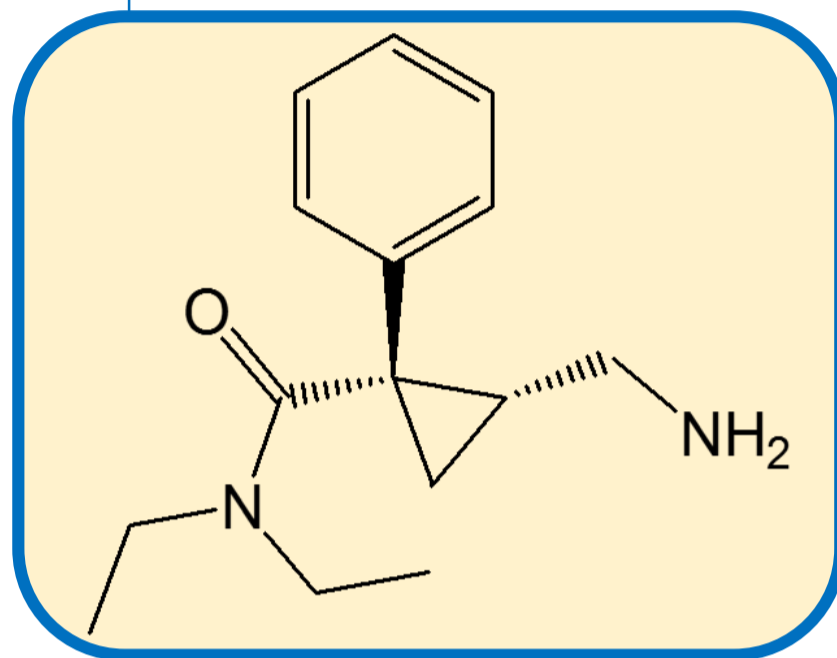

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 9.9	N <sub>1</sub> : 15651	T <sub>1</sub> : 1.3	
tR <sub>2</sub> : 11.5	N <sub>2</sub> : 14830	T <sub>2</sub> : 1.3	Rs: 4.4

