

# DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

## Ruxolitinib

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: Ethanol (100 %)


Flowrate: 0.5 mL/min

Column Oven Temperature: 15°C

UV Detection: 310 nm


Diluent: Ethanol

Sample Concentration: 1.0 mg/mL

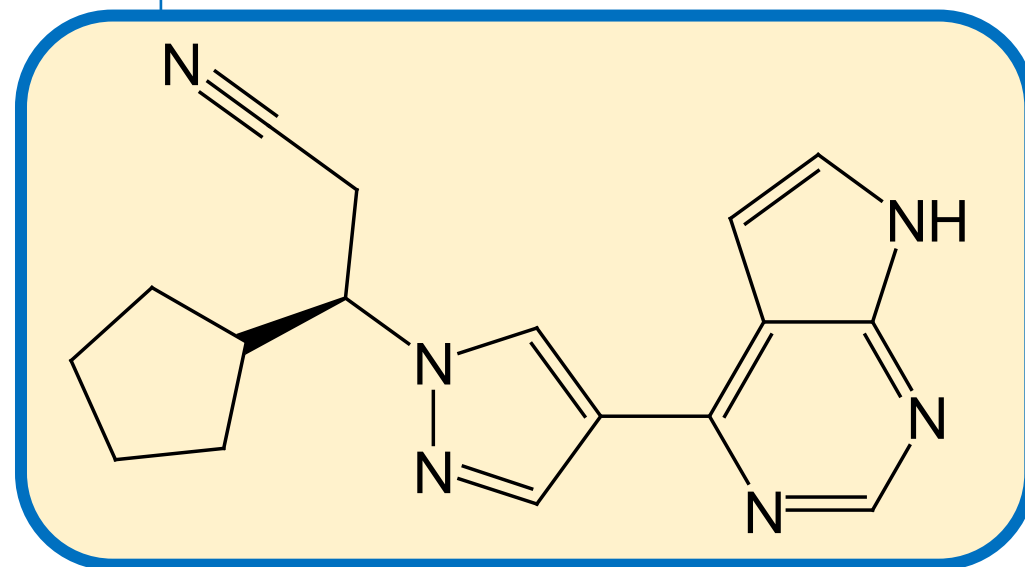

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 16.0	N <sub>1</sub> : 6750	T <sub>1</sub> : 1.2	
tR <sub>2</sub> : 17.9	N <sub>2</sub> : 7800	T <sub>2</sub> : 1.2	Rs: 2.37

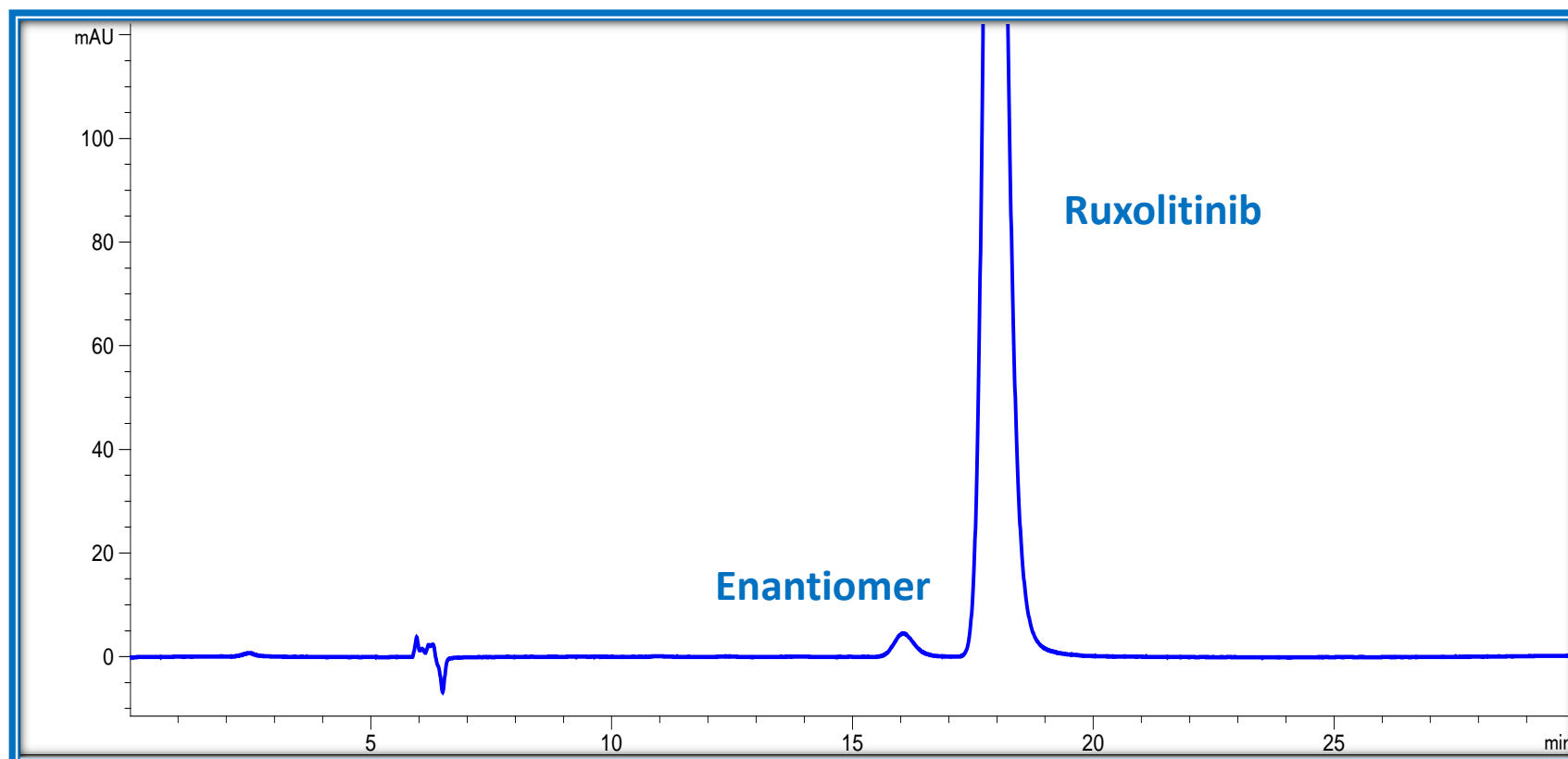

**CHIRAL**  
TECHNOLOGIES  
DAICEL GROUP