

## Tofacitinib

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: Methanol : DEA (100/0.1, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 285 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

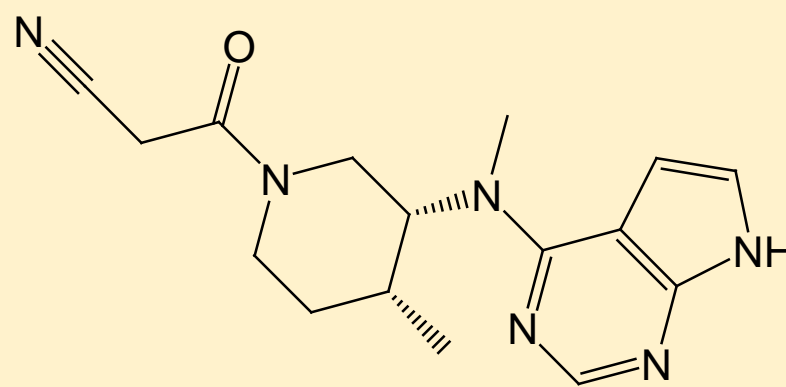

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 11.2	N <sub>1</sub> : 8250	T <sub>1</sub> : 1.0	
tR <sub>2</sub> : 12.6	N <sub>2</sub> : 6200	T <sub>2</sub> : 1.2	Rs: 2.42

