

## Valaciclovir

### CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I(+) [3.0 x 150 mm, 5 µm]

Part No: 53784

Mobile phase: aq. Perchloric acid (pH=1.0) : ACN (80/20, v/v)


Flowrate: 0.2 mL/min

Column Oven Temperature: 25 °C

UV Detection: 254 nm

Diluent: Mobile phase

Sample Concentration: 0.1 mg/mL

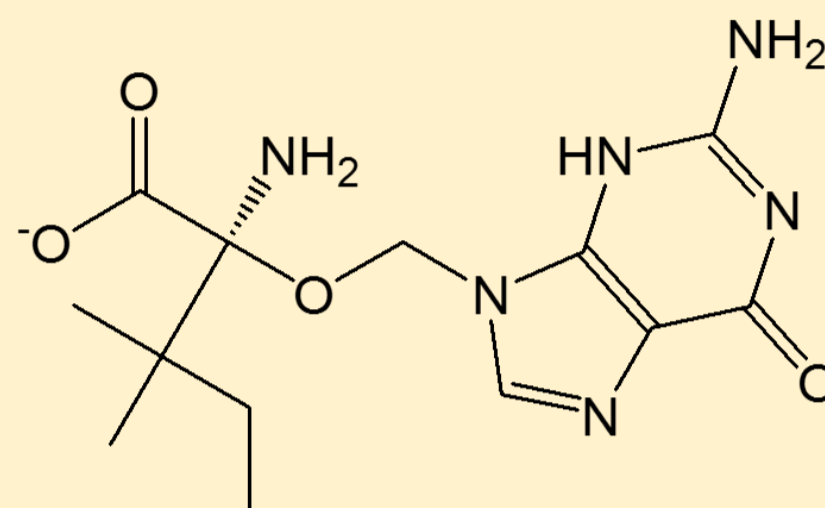

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 5.8	N <sub>1</sub> : 3638	T <sub>1</sub> : 1.5	
tR <sub>2</sub> : 8.2	N <sub>2</sub> : 4756	T <sub>2</sub> : 1.4	Rs: 5.3

