

Valganciclovir

CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I(+) [3.0 x 150 mm, 5 µm]
Part No: 53784

Mobile phase: aq. Perchloric acid (pH=1.0) : ACN (95/05,v/v)


Flowrate: 0.4 mL/min

Column Oven Temperature: 25 °C

UV Detection: 254 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

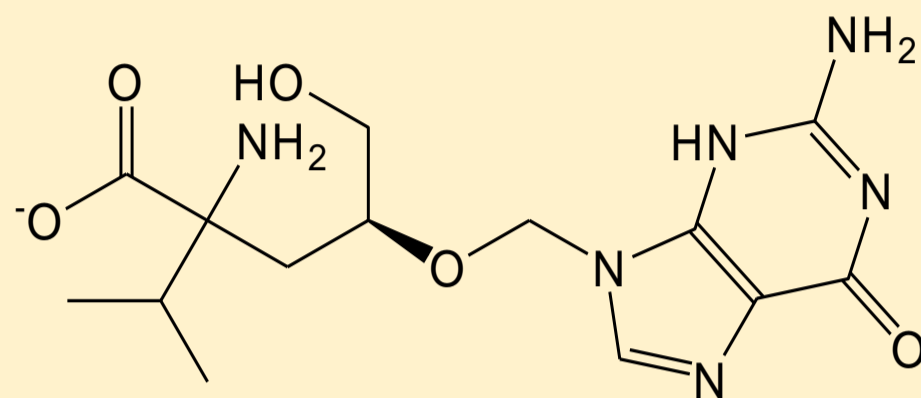

SEPARATION PERFORMANCE

tR ₁ : 3.8	N ₁ : 4040	T ₁ : 1.0	
tR ₂ : 6.2	N ₂ : 2030	T ₂ : 1.6	Rs: 4.0

