

Eliglustat

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IF-3 [4.6 x 250 mm, 3 µm]

Part No: 86525

Mobile phase: n-Hexane : Ethanol : Methanol : DEA
(90/05/05/0.1, v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 280 nm

Diluent: EtOH : n-Hexane (1:1)

Sample Concentration: 10.0 mg/mL

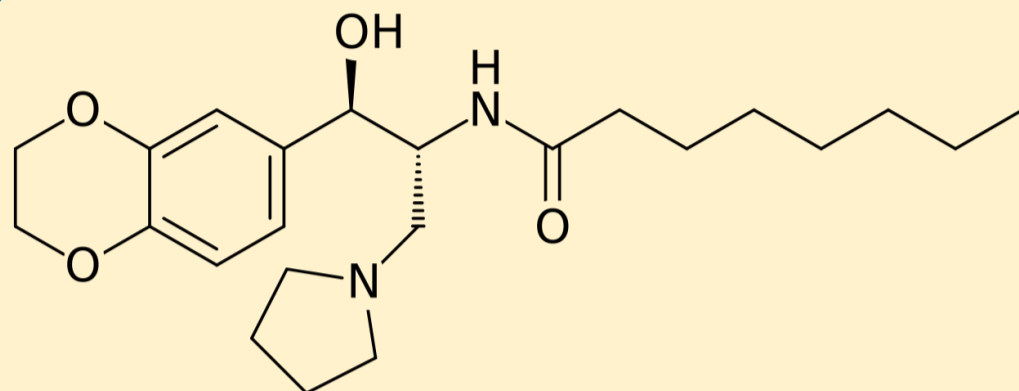

SEPARATION PERFORMANCE

tR ₁ : 13.5	N ₁ : 14290	T ₁ : 1.2	
tR ₂ : 14.9	N ₂ : 13294	T ₂ : 1.2	Rs: 2.5
tR ₃ : 16.0	N ₁ : 13585	T ₁ : 1.2	Rs: 2.1
tR ₄ : 17.8	N ₂ : 5179	T ₂ : 2.3	Rs: 2.3

