

Droxidopa

CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR (+) [4.0 x 150 mm, 5 µm]

Part No: 27714

Mobile phase: aq. Perchloric acid (pH=1.0) (100 %)


Flowrate: 0.2 mL/min

Column Oven Temperature: 10 °C

UV Detection: 280 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

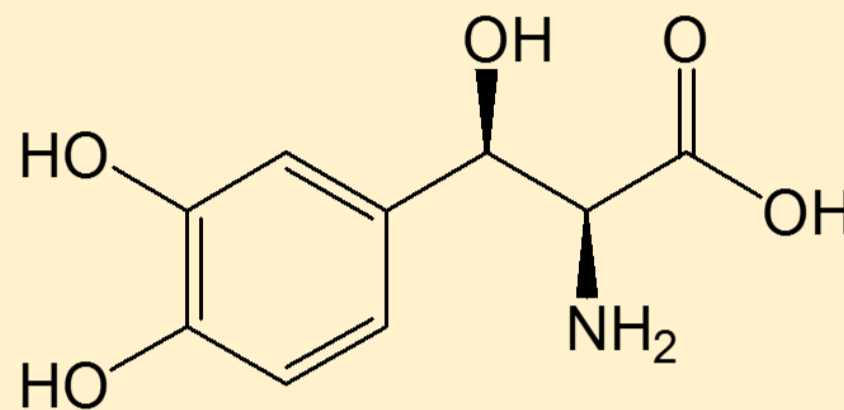

SEPARATION PERFORMANCE

t_{R1} : 10.8	N_1 : 5925	T_1 : 1.3	
t_{R2} : 12.1	N_2 : 5952	T_2 : 1.3	R_s : 2.2

