

Rivaroxaban

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IH [4.6 x 250 mm, 5 µm]

Part No: 89325

Mobile phase: Ethanol : ACN : DEA (90/10/0.1, v/v/v)

Flowrate: 0.7 mL/min


Column Oven Temperature: 40 °C

UV Detection: 250 nm

Sample Concentration: 1.0 mg/mL

Diluent: MeOH : ACN (1 : 1)

Injection Volume : 10 µL

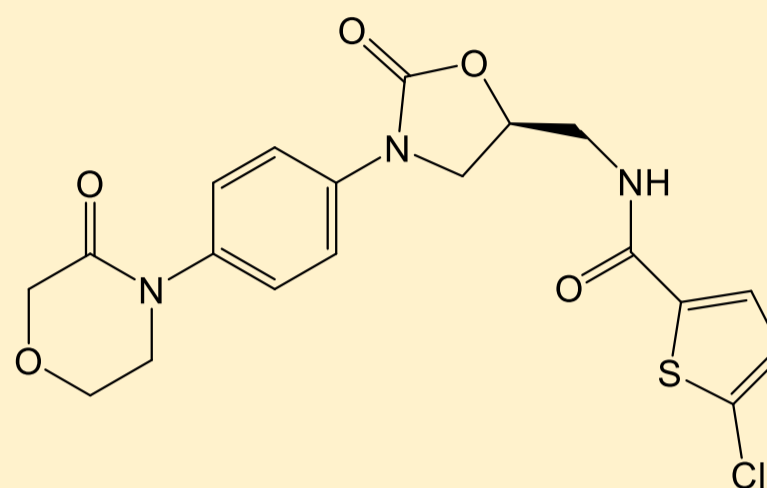

SEPARATION PERFORMANCE

tR_1 : 6.4	N_1 : 9601	T_1 : 1.1	
tR_2 : 7.6	N_2 : 8968	T_2 : 1.2	R_s : 4.2

