

Brinzolamide

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 150 mm, 3 µm]

Part No: 87524

Mobile phase: aq. 10 mM Ammonium bicarbonate : ACN (70/30, v/v)

Flowrate: 1.0 mL/min


Column Oven Temperature: 25 °C

UV Detection: 254 nm

Sample Concentration: 1.0 mg/mL

Diluent: Mobile phase

Injection Volume: 10 µL

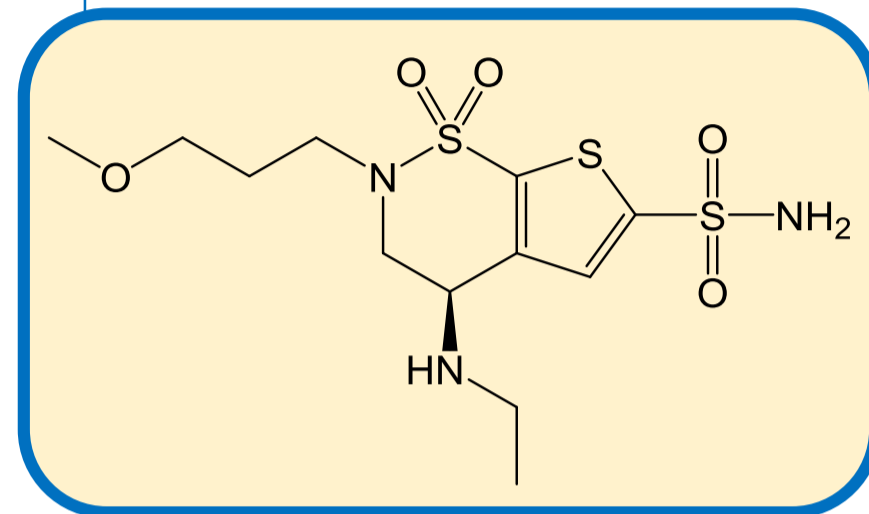

SEPARATION PERFORMANCE

t_{R1} : 5.10

N_1 : 5086

T_1 : 1.1

t_{R2} : 6.23

N_2 : 2865

T_2 : 1.1

R_s : 2.99

