

Baclofen

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ZWIX(+) [3.0 x 150mm, 3µm]

Part No: 51584

Mobile phase: 5 mM Formic Acid + 5 mM ammonium formate in
(Methanol:ACN:H₂O 75/20/05, v/v/v)


Flowrate: 0.7 mL/min

Column Oven Temperature: 25 °C

UV detection: 230 nm

Diluent: H₂O : Methanol (2:8)

Concentration: 0.7 mg/mL.

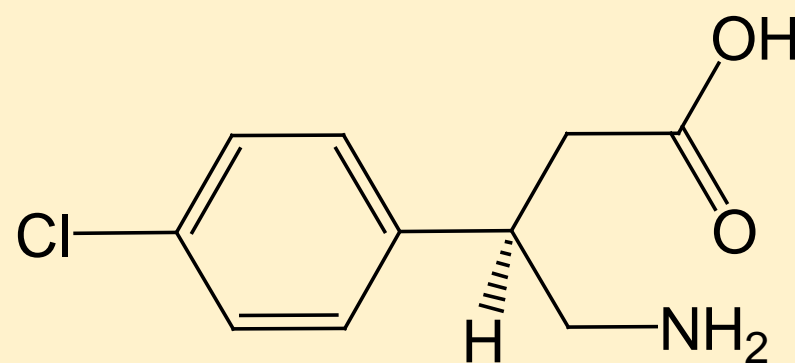

SEPARATION PERFORMANCE

tR ₁ : 4.44	N ₁ : 5900	T ₁ : 1.5	
tR ₂ : 5.45	N ₂ : 6100	T ₂ : 1.5	Rs: 4.02

