

Baclofen

CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I(+) [3.0 x 150 mm, 5 µm]

Part No: 53784

Mobile phase: aq. Perchloric acid (pH=1.5):ACN (60/40,v/v)


Flowrate: 0.4 mL/min

Column Oven Temperature: 25 °C

UV Detection: 230 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

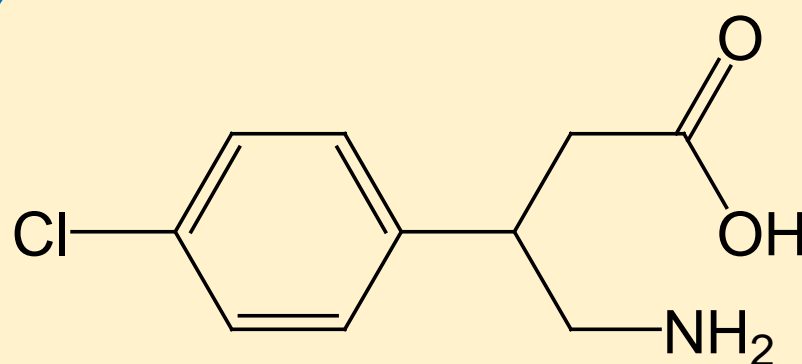

SEPARATION PERFORMANCE

tR ₁ : 3.1	N ₁ : 3964	T ₁ : 1.1	
tR ₂ : 4.9	N ₂ : 4402	T ₂ : 1.4	Rs: 7.15

