

Batefenterol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ZWIX(-) [3.0 x 250 mm, 3 μm]
Part No: 52585

Mobile phase: 50 mM formic acid + 50 mM DEA in
(Methanol : THF : H₂O 50/50/01, v/v/v)


Flowrate: 0.5 mL/min

Column Oven Temperature: 25 °C

UV Detection: 254 nm

Diluent: Methanol

Sample Concentration: 1.0 mg/mL

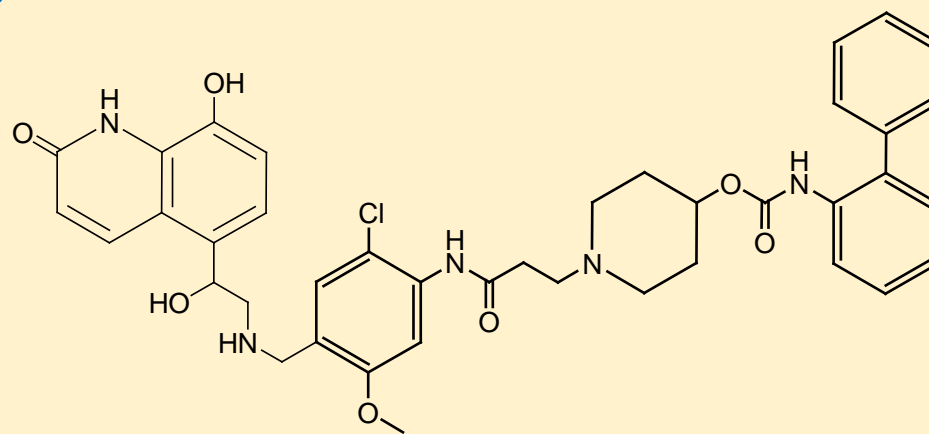

SEPARATION PERFORMANCE

tR ₁ : 4.0	N ₁ : 5400	T ₁ : 1.9	
tR ₂ : 5.3	N ₂ : 4671	T ₂ : 2.0	Rs: 5.0

