


Canagliflozin

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 150 mm, 5 µm]
 Part No: 87324
 Mobile phase: M.P. A: 10 mM KPF6 (pH=2.0), M.P. B: ACN
 T/%B: 0/30, 15/90, 15.1/30, 20/30
 Flowrate: 1.0 mL/min
 Column Oven Temperature: 25 °C
 UV Detection: 290 nm
 Diluent: Mobile phase
 Sample Concentration: 1.0 mg/mL

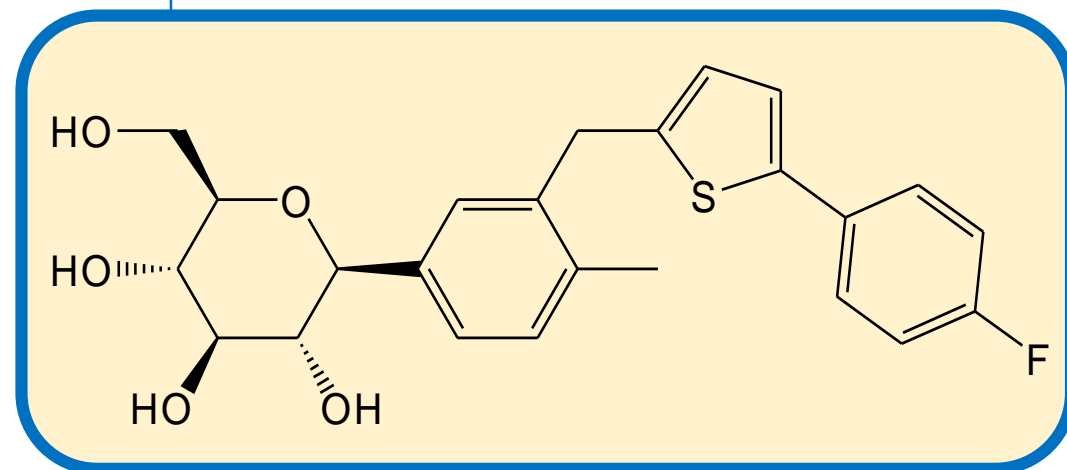

SEPARATION PERFORMANCE

tR ₁ : 8.1	N ₁ : 12429	T ₁ : 1.2	
tR ₂ : 9.8	N ₂ : 14511	T ₂ : 1.2	Rs: 5.7
tR ₃ : 14.6	N ₃ : 7018	T ₃ : 1.1	Rs: 9.4

