

## Cinacalcet

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® AD-3R [4.6 x 150 mm, 3 µm]

Part No: 19824

Mobile phase: 5mM ammonium bicarbonate in  
(H<sub>2</sub>O:Methanol 10/90, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 225 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

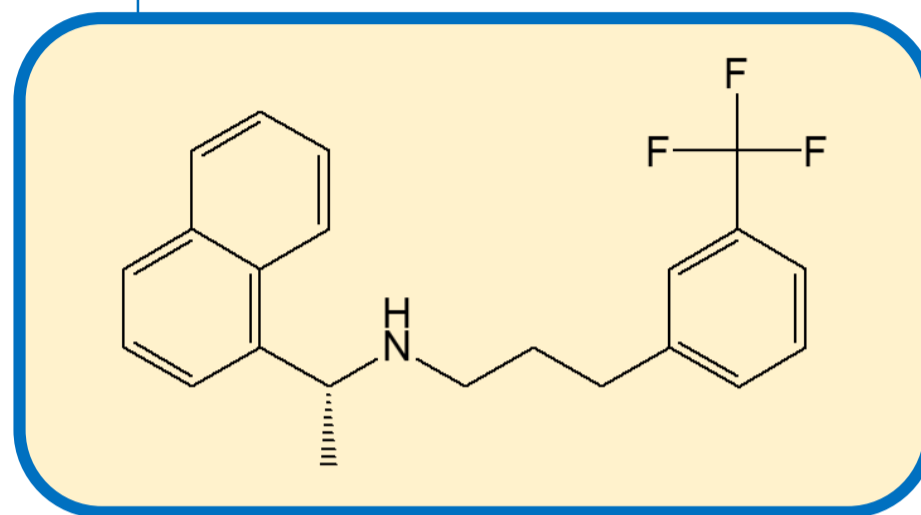

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 3.6	N <sub>1</sub> : 5052	T <sub>1</sub> : 1.3	
tR <sub>2</sub> : 4.7	N <sub>2</sub> : 4800	T <sub>2</sub> : 1.4	Rs: 2.9

