

Dexrazoxane

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK®ID [4.6 x 250 mm, 5 µm]

Part No: 84325

Mobile phase: Methanol : ACN : DEA (60:40:0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV detection: 235 nm

Diluent: Mobile phase

Concentration: 1.0 mg/mL.

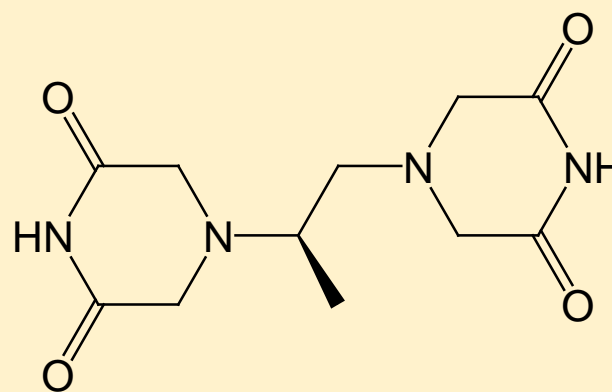

SEPARATION PERFORMANCE

| | | | |
|------------------------|-----------------------|----------------------|----------|
| tR ₁ : 3.68 | N ₁ : 9500 | T ₁ : 1.4 | |
| tR ₂ : 5.00 | N ₂ : 6100 | T ₂ : 2.4 | Rs: 6.47 |

