

Dropropizine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IF-3 [4.6 x 250 mm, 3 µm]

Part No: 86525

Mobile phase: n-Hexane : EtOH : EDA (80/20/0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 250 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

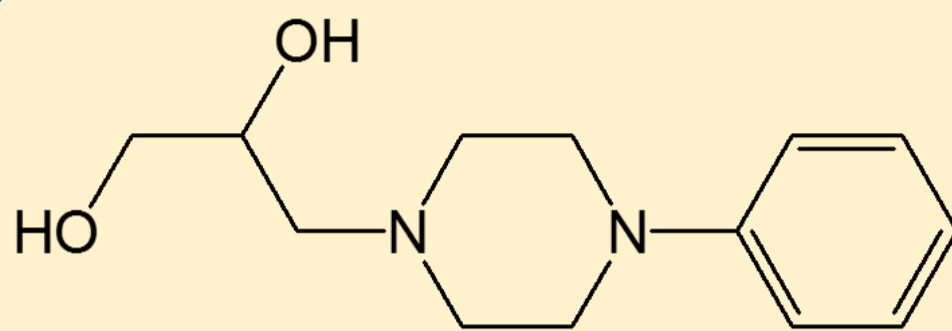

SEPARATION PERFORMANCE

tR_1 : 13.7	N_1 : 22030	T_1 : 1.4	
tR_2 : 15.5	N_2 : 21806	T_2 : 1.4	R_s : 4.41

