

Flupentixol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ID-3 [4.6 x 250 mm, 3 µm]
Part No: 84525

Mobile phase: n-Hexane : 2-Butanol : Methanol : DEA
(93/07/0.5/0.3, v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 15 °C

UV Detection: 260 nm

Diluent: Methanol : 2-propanol : n-Hexane (2:2:6)

Sample Concentration: 0.5 mg/mL

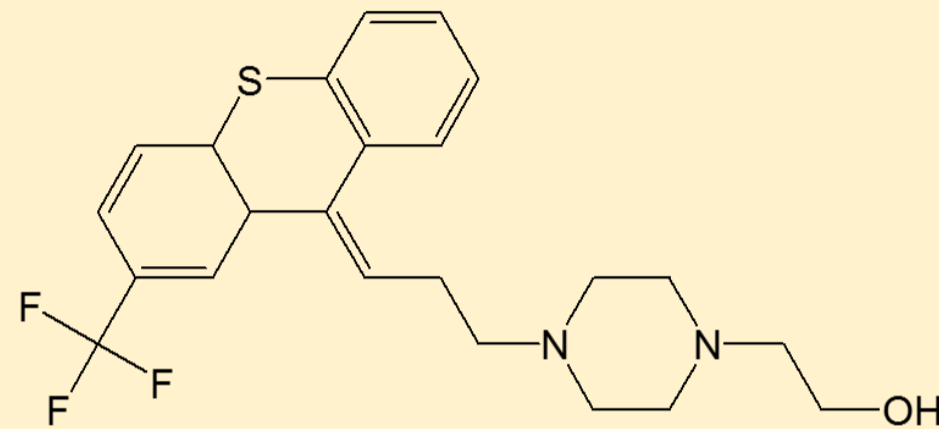

SEPARATION PERFORMANCE

t_{R1} : 12.9	N_1 : 7098	T_1 : 1.6	
t_{R2} : 15.8	N_2 : 2909	T_2 : 1.4	R_s : 3.3

