

Ivabradine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 250 mm, 5 µm]

Part No: 87325

Mobile phase: Methanol : ACN : DEA (50/50/0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 285 nm

Diluent: Mobile phase

Sample Concentration: 2.0 mg/mL

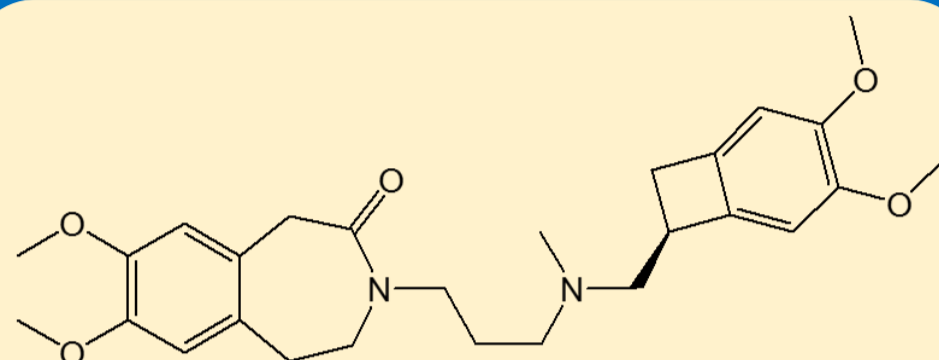

SEPARATION PERFORMANCE

t_{R1} : 9.4	N_1 : 5800	T_1 : 1.0	
t_{R2} : 12.1	N_2 : 3100	T_2 : 1.6	R_s : 4.0

