

## Linezolid

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IA-3 [4.6 x 250 mm, 3 µm]

Part No: 80525

Mobile phase: ACN : Methanol : DEA (90:10:0.1, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 260 nm

Diluent: Mobile phase

Sample Concentration: 0.25 mg/mL

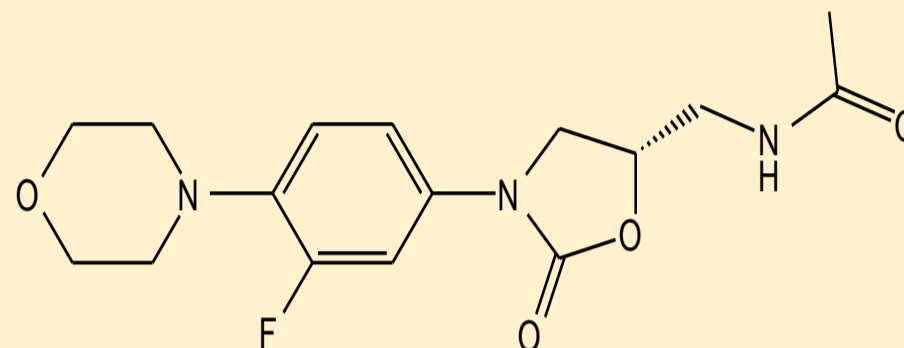

### SEPARATION PERFORMANCE

$t_{R1}$ : 4.1	$N_1$ : 6806	$T_1$ : 2.2	
$t_{R2}$ : 5.0	$N_2$ : 3911	$T_2$ : 1.9	$R_s$ : 3.8

