

Linezolid

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IA-3 [4.6 x 150 mm, 3 µm]
Part No: 80524

Mobile phase: 5mM ammonium bicarbonate in
(H₂O : ACN, 05/95, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 260 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

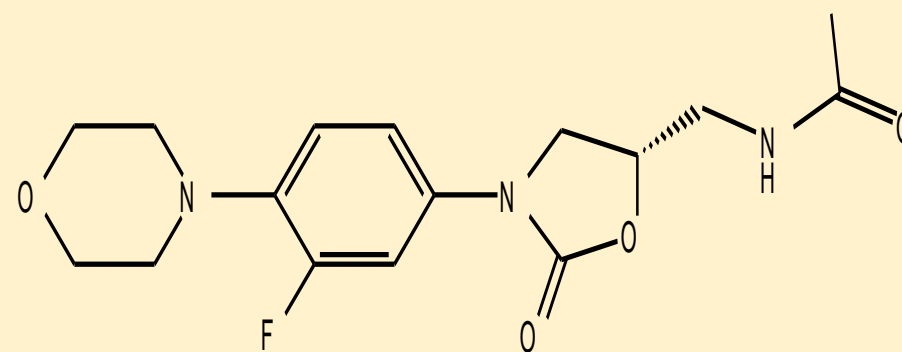

SEPARATION PERFORMANCE

tR ₁ : 2.5	N ₁ : 7535	T ₁ : 1.4	
tR ₂ : 3.1	N ₂ : 3719	T ₂ : 1.6	Rs: 3.95

