

Lisdexamfetamine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IF-3 [4.6 x 250 mm, 3 µm]
Part No: 86525

Mobile phase: n-Hexane : MTBE : Methanol : EDA
(30/30/40/0.1,v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 258 nm

Diluent: Methanol : MTBE : n-Hexane (4:3:3)

Sample Concentration: 2.0 mg/mL

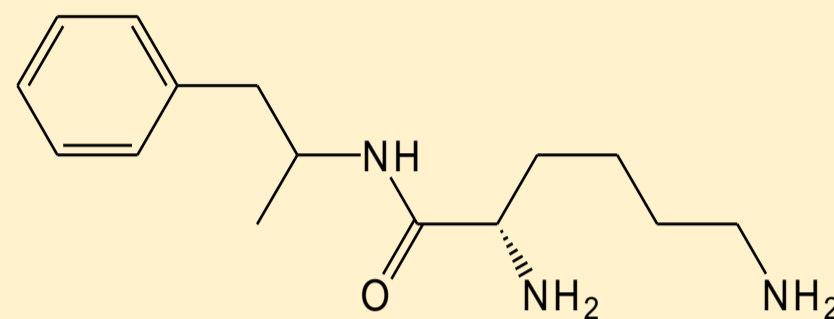

SEPARATION PERFORMANCE

t_{R1} : 5.2	N_1 : 5582	T_1 : 1.9	
t_{R2} : 6.8	N_2 : 9557	T_2 : 1.5	R_s : 5.8

