

Mirabegron

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IF-3 [4.6 x 250 mm, 3 µm]
Part No: 86525

Mobile phase: n-Hexane : MTBE : Methanol : EDA
(30/30/40/0.1, v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 250 nm

Diluent: Mobile phase

Sample Concentration: 0.2 mg/mL

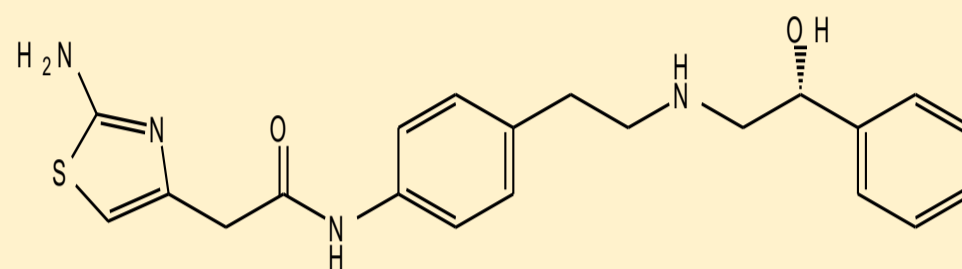

SEPARATION PERFORMANCE

t _{R1} : 8.7	N ₁ : 6893	T ₁ : 1.3	
t _{R2} : 10.1	N ₂ : 8066	T ₂ : 1.4	Rs: 3.1

