

Mitiglinide

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IE [4.6 x 250 mm, 5 µm]

Part No: 85325

Mobile phase: MtBE : ACN : TFA (80/20/0.3, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 215 nm

Diluent: Methanol : Mobile phase (1 : 9)

Sample Concentration: 1.0 mg/mL

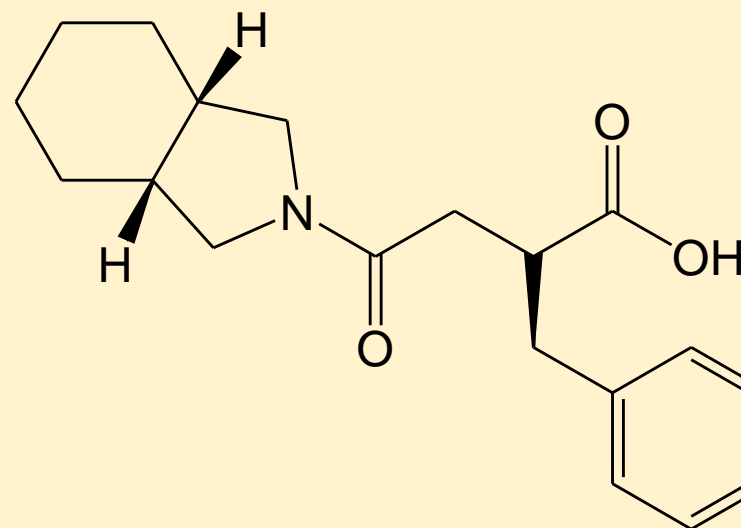

SEPARATION PERFORMANCE

t_{R1} : 8.9	N_1 : 4026	T_1 : 1.2	
t_{R2} : 10.2	N_2 : 6201	T_2 : 1.2	R_s : 2.4

