


Ranolazine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 150 mm, 5 µm]
Part No: 87324
Mobile phase: 10 mM ammonium bicarbonate in
(H₂O:Methanol (5/95, v/v))
Flowrate: 1.0 mL/min
Column Oven Temperature: 40 °C
UV Detection: 220 nm
Diluent: Mobile phase
Sample Concentration: 1.0 mg/mL

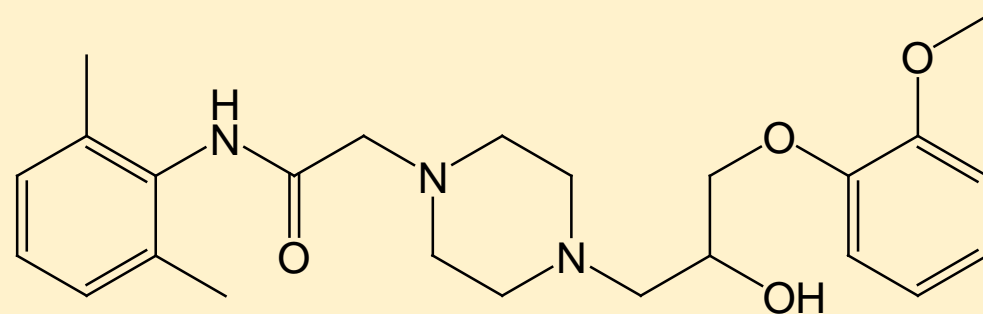

SEPARATION PERFORMANCE

tR ₁ : 5.19	N ₁ : 3909	T ₁ : 1.0	
tR ₂ : 6.90	N ₂ : 3136	T ₂ : 1.1	Rs: 4.12

