

Salmeterol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ZWIX(+) [3.0 x 150 mm, 3 µm]

Part No: 51584

Mobile phase: 12.5 mM 2-Amino-1-butanol in
(Methanol : ACN: H₂O 49/49/02, v/v/v)


Flowrate: 0.5 mL/min

Column Oven Temperature: 25 °C

UV Detection: 225 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

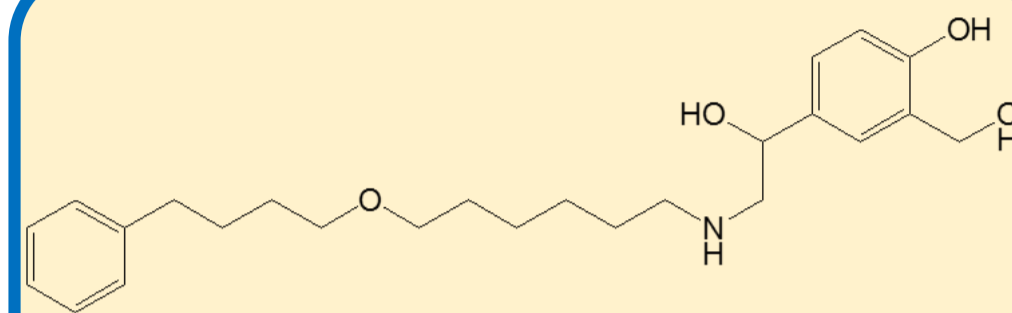

SEPARATION PERFORMANCE

tR ₁ : 9.8	N ₁ : 6485	T ₁ : 1.3	
tR ₂ : 11.6	N ₂ : 6832	T ₂ : 1.4	Rs: 3.43

