

Thalidomide

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® AD-3R [4.6 x 150 mm, 3 µm]
Part No: 19824

Mobile phase: 5 mM ammonium acetate in
(H₂O : Methanol : ACN 05/25/70, v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 220 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

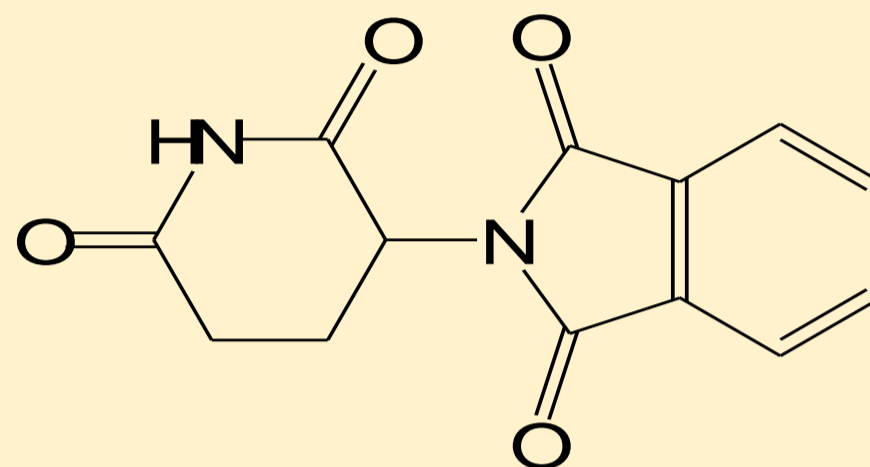

SEPARATION PERFORMANCE

t _{R1} : 3.5	N ₁ : 6640	T ₁ : 1.2	
t _{R2} : 4.3	N ₂ : 4030	T ₂ : 1.3	Rs: 3.3

