

Zopiclone

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IC-3 [4.6 x 150 mm, 5 µm]
Part No: 83524

Mobile phase: 5mM ammonium bicarbonate in
(H₂O:ACN 05/95, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 304 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

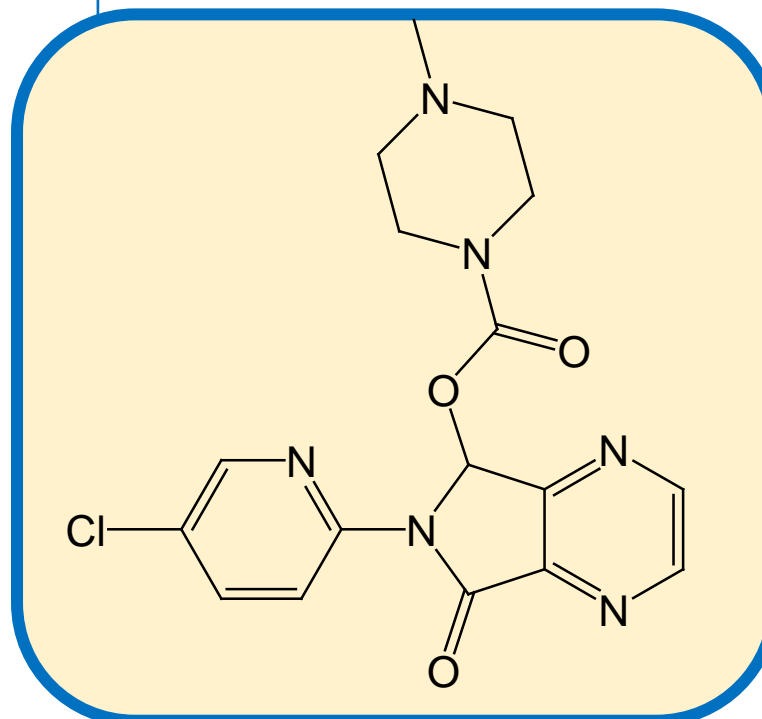

SEPARATION PERFORMANCE

t _{R1} : 3.7	N ₁ : 10438	T ₁ : 1.0	
t _{R2} : 5.0	N ₂ : 10555	T ₂ : 1.1	Rs: 8.0

