

Atazanavir

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® ID-3 [4.6 x 250 mm, 3 µm]

Part No: 84525

Mobile phase: Methanol : DEA (100/0.1, v/v)


Flowrate: 0.7 mL/min

Column Oven Temperature: 25 °C

UV Detection: 280 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

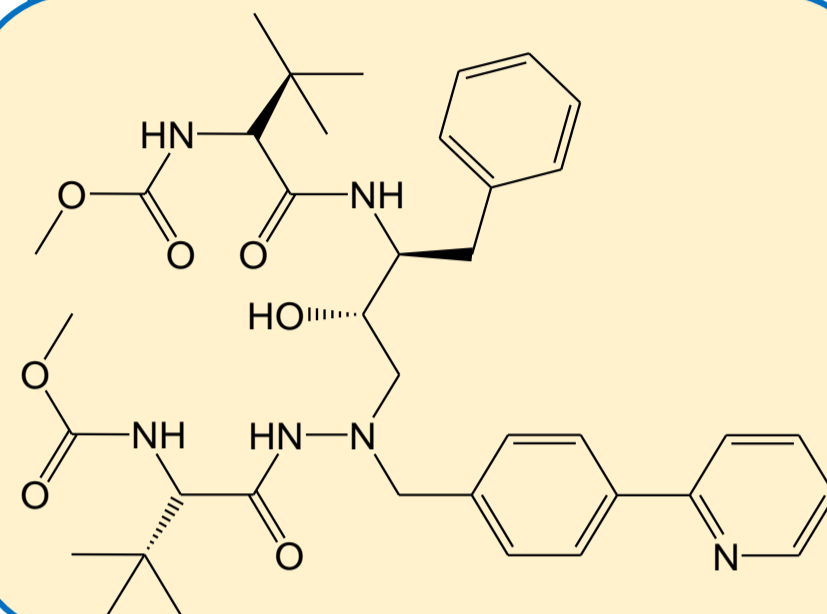

SEPARATION PERFORMANCE

tR_1 : 6.6	N_1 : 4206	T_1 : 1.8	
tR_2 : 8.2	N_2 : 3872	T_2 : 2.0	R_s : 3.3

