

## Mosapride

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK®IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: n-Hexane : Ethanol : Methanol : DEA (80/10/10/0.1, v/v/v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 270 nm

Diluent: EtOH

Sample Concentration: 1.0 mg/mL

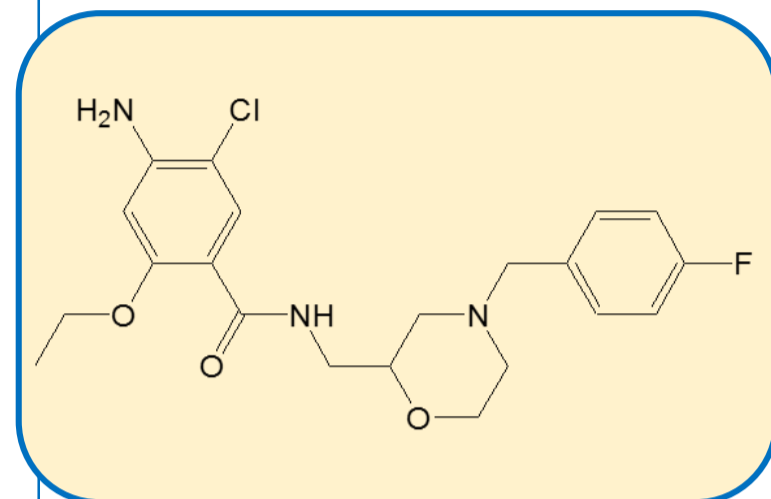

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 17.1	N <sub>1</sub> : 9145	T <sub>1</sub> : 1.7	
tR <sub>2</sub> : 18.8	N <sub>2</sub> : 9554	T <sub>2</sub> : 1.7	Rs: 2.3

