

Pidotimod

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: Methanol : ACN : TFA (60/40/0.1, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 230 nm

Diluent: Mobile phase

Sample Concentration: 2.0 mg/mL

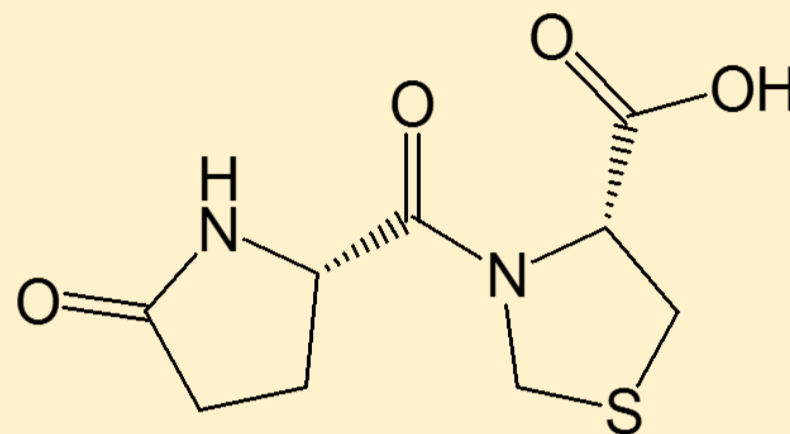

SEPARATION PERFORMANCE

t_{R1} : 5.1	N_1 : 9038	T_1 : 1.7	
t_{R2} : 6.2	N_2 : 3474	T_2 : 1.4	R_s : 3.6

