

Propranolol

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 150 mm, 3 µm]

Part No: 87524

Mobile phase: 20 mM ammonium bicarbonate in
(H₂O : MeOH (20/80, v/v))


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 290 nm

Diluent: Mobile phase

Sample Concentration: 0.1 mg/mL

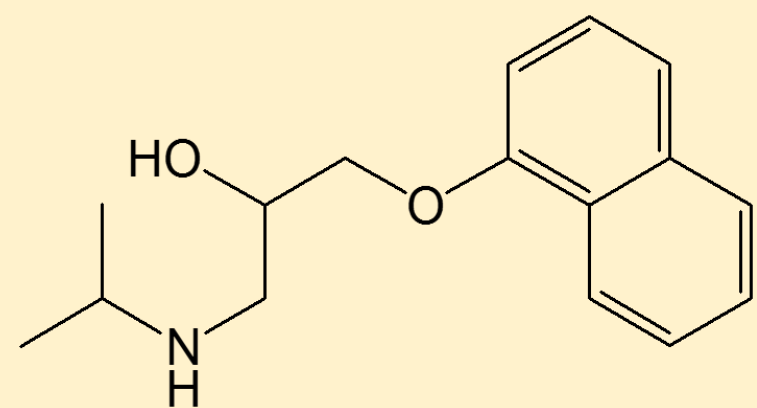

SEPARATION PERFORMANCE

tR_1 : 3.6	N_1 : 4700	T_1 : 1.4	
tR_2 : 4.4	N_2 : 4215	T_2 : 1.2	R_s : 3.6

