

Rabeprazole

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IB-3 [4.6 x 150 mm, 3 µm]

Part No: 81524

Mobile phase: 5 mM ammonium acetate in (H₂O : ACN 5/95, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 280 nm

Diluent: Mobile phase

Sample Concentration: 0.5 mg/mL

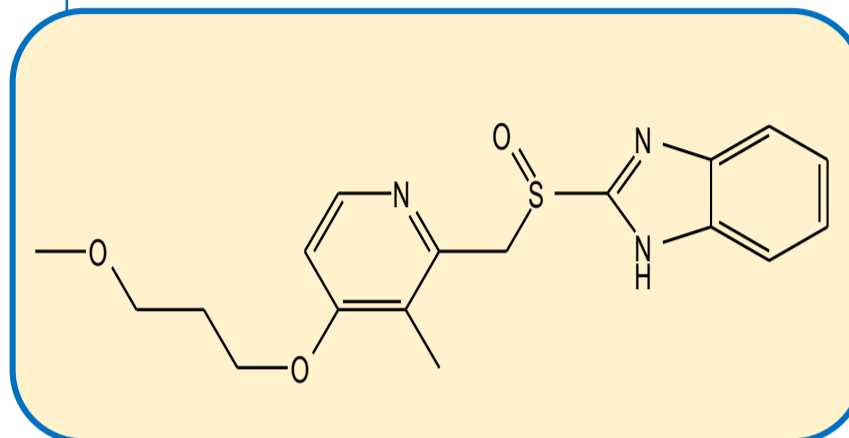

SEPARATION PERFORMANCE

tR ₁ : 3.1	N ₁ : 3354	T ₁ : 1.3	
tR ₂ : 4.0	N ₂ : 2082	T ₂ : 1.3	Rs: 3.0

