

## Voriconazole

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IC-3 [4.6 x 150 mm, 3 µm]

Part No: 83524

Mobile phase: aq. 10mM ammonium bicarbonate : ACN (50/50, v/v)


Flowrate: 1.0 mL/min

Column Oven Temperature: 25 °C

UV Detection: 254 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

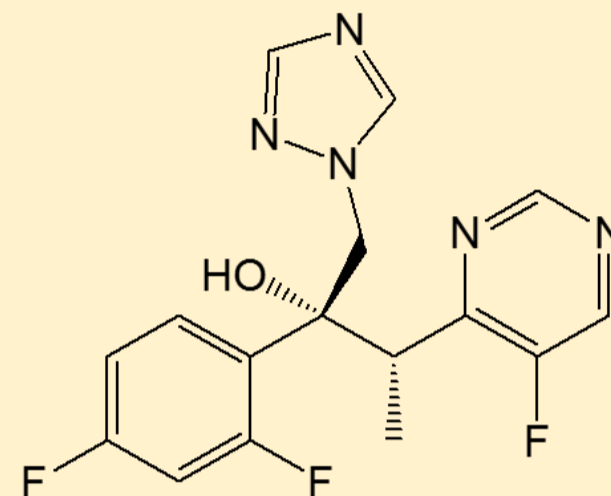

### SEPARATION PERFORMANCE

$tR_1$ : 9.0	$N_1$ : 4970	$T_1$ : 1.1	
$tR_2$ : 11.0	$N_2$ : 3844	$T_2$ : 1.4	$R_s$ : 3.2

