

Mirabegron

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG [4.6 x 150 mm, 5 µm]

Part No: 87324

Mobile phase: 5 mM ammonium acetate in
(H₂O : Methanol (5/95, v/v))


Flowrate: 1.0 mL/min

Column Oven Temperature: 40 °C

UV Detection: 250 nm

Diluent: Mobile phase

Sample Concentration: 1.0 mg/mL

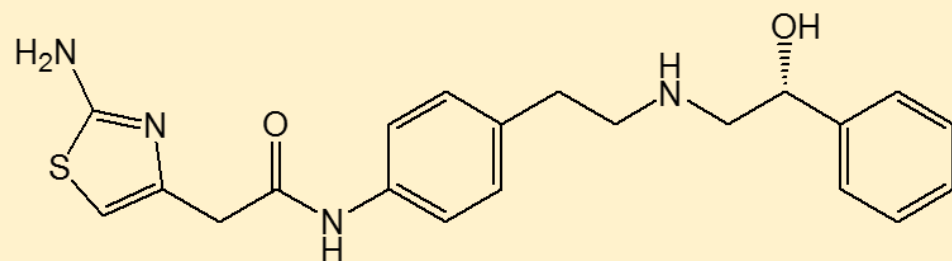

SEPARATION PERFORMANCE

tR ₁ : 5.2	N ₁ : 3759	T ₁ : 1.1	
tR ₂ : 6.5	N ₂ : 3459	T ₂ : 1.1	Rs: 3.24

