

## Safinamide

### CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IH-3 [4.6 x 150 mm, 3 µm]  
Part No: 89524

Mobile phase: 5 mM Ammonium bicarbonate in  
[(Water : ACN) (05/95, v/v)]

Flowrate: 0.7 mL/min


Column Oven Temperature: 15 °C

UV Detection: 226 nm

Sample Concentration: 1.0 mg/mL

Diluent: Mobile phase

Injection Volume: 5 µL

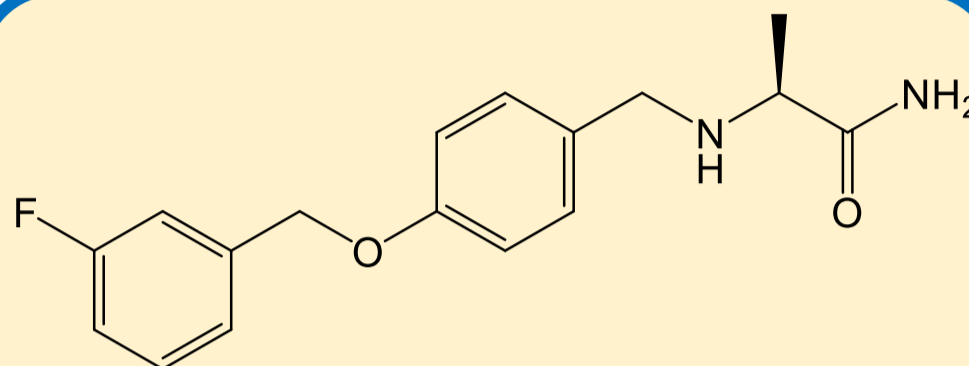

### SEPARATION PERFORMANCE

$t_{R1}$ : 4.3	$N_1$ : 15458	$T_1$ : 1.3	
$t_{R2}$ : 4.9	$N_2$ : 15877	$T_2$ : 1.2	$R_s$ : 4.5

