

Frovatriptan

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 150 mm, 3 µm]
Part No: 87524

Mobile phase: 10 mM ammonium bicarbonate in
[(Water: ACN) (10/90, v/v)]

Flowrate: 1.0 mL/min


Column Oven Temperature: 40 °C

UV Detection: 245 nm

Sample Concentration: 0.2 mg/mL

Diluent: Mobile phase

Injection Volume: 10 µL

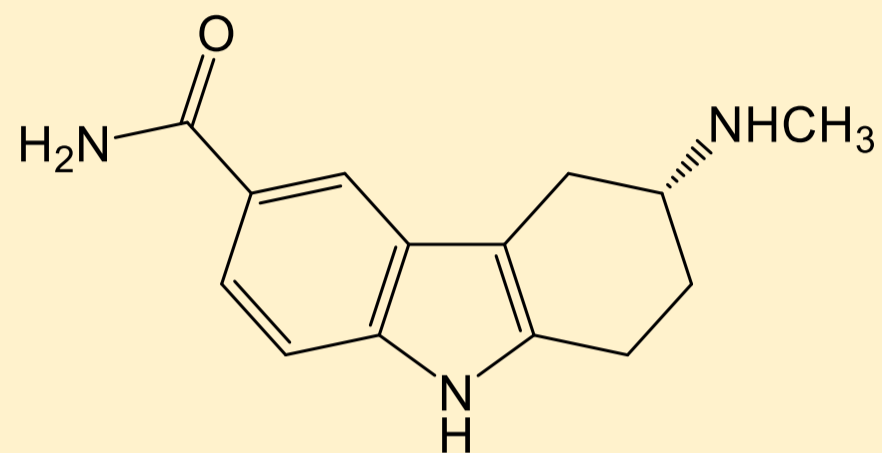

SEPARATION PERFORMANCE

t_{R_1} : 4.8	N_1 : 11751	T_1 : 1.4	
t_{R_2} : 5.3	N_2 : 8464	T_2 : 1.9	R_s : 2.3

