

Verapamil

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 250 mm, 3 µm]

Part No: 87525

Mobile phase: n-Hexane : 2-propanol : DEA (70/30/0.1, v/v/v)

Flowrate: 1.0 mL/min


Column Oven Temperature: 25 °C

UV Detection: 230 nm

Sample Concentration: 2.0 mg/mL

Diluent: Ethanol : n-Hexane (1 : 1)

Injection Volume: 10µL

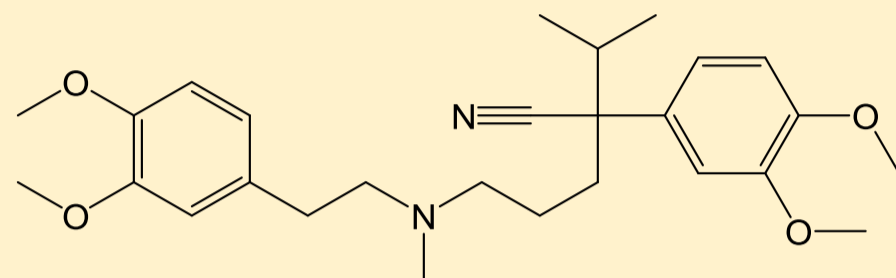

SEPARATION PERFORMANCE

t_{R1} : 11.0	N_1 : 7015	T_1 : 1.5	
t_{R2} : 12.6	N_2 : 6666	T_2 : 1.5	R_s : 2.7

