

Prilocaine

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IG-3 [4.6 x 150 mm, 3 µm]

Part No: 87524

Mobile phase: 5 mM Ammonium bicarbonate in
[(Water : ACN) (5/95, v/v)]

Flowrate: 0.7 mL/min


Column Oven Temperature: 25 °C

UV Detection: 240 nm

Sample Concentration: 1.0 mg/mL

Diluent: Mobile phase

Injection Volume: 10 µL

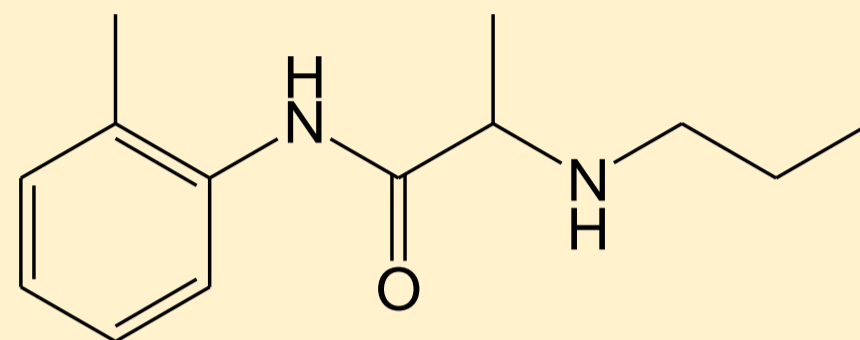

SEPARATION PERFORMANCE

t_{R1} : 4.3	N_1 : 11910	T_1 : 1.2	
t_{R2} : 5.5	N_2 : 9393	T_2 : 1.2	R_s : 6.1

