

## Indoximod

### CHROMATOGRAPHIC CONDITIONS

Column: CROWNPAK® CR-I (+) (3.0 x 150 mm, 5 µm)  
Part No: 53784

Mobile phase: aq. HClO<sub>4</sub> ( pH: 2.0 ) : ACN ( 80/20, v/v)

Flow rate: 0.4 mL/min


Column Oven Temperature: 25°C

UV Detection: 210 nm

Sample Concentration: 0.5 mg/mL

Diluent: aq. HClO<sub>4</sub> ( pH: 1.0 )

Injection Volume: 10µL

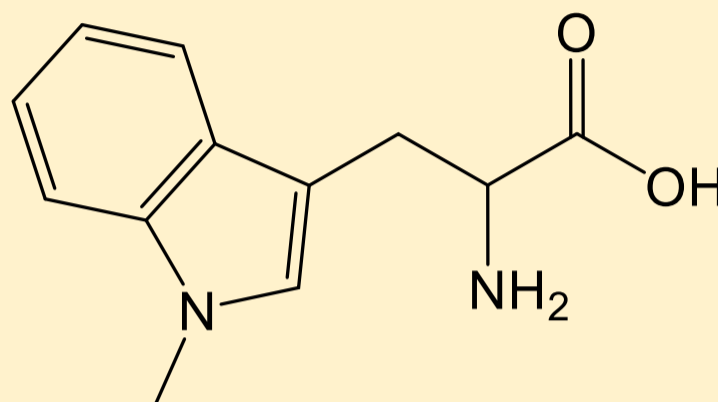

### SEPARATION PERFORMANCE

tR <sub>1</sub> : 2.8	N <sub>1</sub> : 4585	T <sub>1</sub> : 1.4			
tR <sub>2</sub> : 4.2	N <sub>2</sub> : 4517	T <sub>2</sub> : 1.5	α: 2.2	Rs: 6.8	

