

REMDESIVIR Pentafluoro Phosphoramidite Key Starting Material

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK® IA-3 [250 mm Xx 4.6 mm i.d., 3.0µm]

Part No. 80525

Mobile phase : n-Hexane/IPA = 90/10 (v/v)

Flowrate: 1.0 ml/min


Column Oven Temperature: 25°C

UV Detection: 210 nm

Sample Concentration: 1.0 mg/ml

Diluent: n-Hexane/IPA = 90/10 (v/v)

Injection Volume: 10 µl

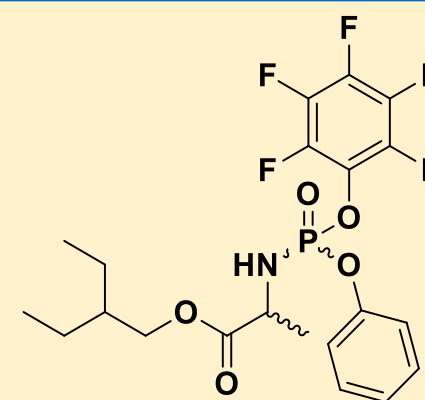

Remdesivir Pentafluoro Phosphoramidite
Key Starting Material

tR ₁ : 7.0	N ₁ : 20048	T ₁ : 1.1	
tR ₂ : 7.6	N ₂ : 21219	T ₂ : 1.3	Rs: 2.8
tR ₃ : 8.1	N ₂ : 26014	T ₂ : 1.2	Rs: 2.6
tR ₄ : 11.6	N ₂ : 18338	T ₂ : 1.3	Rs: 12.8

