

## DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

### Baloxavir

Column: CHIRALPAK® IH-3 [4.6 x 250 mm, 3 µm]

Part No: 89525

Mobile Phase: Methyl tert-butyl ether : Ethanol: Ethylenediamine (50/50/0.1,v/v/v)


Column Oven Temperature: 40°C

UV Detection: 260nm

Flow Rate: 1.0 ml/min

Diluent: Mobile Phase

Sample Concentration: 1.0 mg/ml

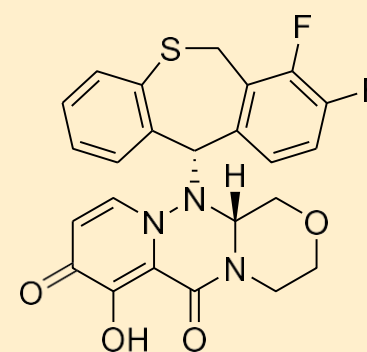

tR <sub>1</sub> : 5.8	N <sub>1</sub> : 50646	T <sub>1</sub> :1.2	
tR <sub>2</sub> : 7.4	N <sub>2</sub> : 45127	T <sub>2</sub> :1.0	Rs: 4.9
tR <sub>3</sub> : 7.9	N <sub>3</sub> : 45134	T <sub>3</sub> :1.2	Rs: 1.3
tR <sub>4</sub> : 12.0	N <sub>4</sub> : 17202	T <sub>4</sub> :1.3	Rs: 6.3

