

DAICEL CHIRAL TECHNOLOGIES GLOBAL APPLICATION CENTRE

Sofosbuvir

Column: CHIRALPAK® ID-3 [4.6 x 250 mm, 3 µm]

Part No: 84525

Mobile Phase: n-Hexane : 2-Propanol : 2-Butanol : DEA (50/25/25/0.1,v/v/v/v)


Column Oven Temperature: 40°C

UV Detection: 260nm

Flow Rate: 1.0 ml/min

Diluent: Mobile Phase

Sample Concentration: 1.0 mg/ml

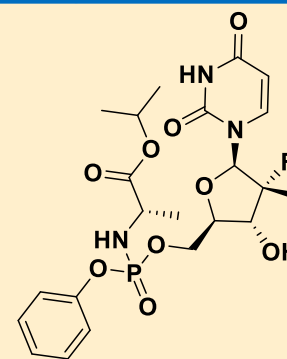

tR ₁ : 8.4	N ₁ : 5190	T ₁ :1.7	
tR ₂ : 10.5	N ₂ : 5316	T ₂ :1.2	Rs: 4.2
tR ₃ : 13.0	N ₃ : 5323	T ₃ :1.3	Rs: 3.8
tR ₄ : 15.2	N ₄ : 3942	T ₄ :1.3	Rs: 2.6

