


Flurbiprofen

CHROMATOGRAPHIC CONDITIONS

Column: CHIRALPAK®IG [250 mm x 4.6 mm i.d., 5.0 μm]
(Part No. 87325)
Mobile phase: n-Hexane/ Ethanol/Methanol (80/10/10 v/v/v)
Flowrate: 1.0 mL/min
Column Oven Temperature: 25 °C
UV Detection: 220 nm
Diluent: Mobile Phase
Sample Concentration: 0.5 mg/mL
Injection Volume: 10 μL

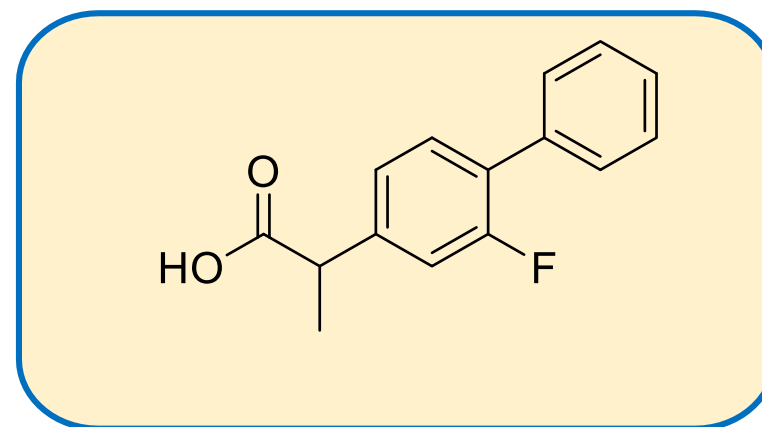

SEPARATION PERFORMANCE

t_{R1} : 4.8	N_1 : 4566	T_1 : 1.4	
t_{R2} : 6.3	N_2 : 5017	T_2 : 1.4	R_s : 4.6

